

Lab Notebook for

PHY 445

Will Lancer

`will.m.lancer@gmail.com`

Lab notebook for PHY 445 Spring 2026

Contents

1 Feb 5, 2026	3
1.1 Experiment: Nuclear Magnetic Resonance	3
2 Feb 10, 2026	4
2.1 Sub-experiment: mineral oil	4
2.2 Enter Professor Du, the King of NMR™ 	4
3 Feb 12, 2026	5
3.1 Beeswax T_1 measurement	5
3.2 Nevermind, fml	5
3.3 Beeswax, pt. 2	6
4 Feb 17, 2026	8
4.1 Spin-echo: light mineral oil	8
4.2 Beeswax-maxxing	9
5 Feb 19, 2026	10
5.1 Light mineral oil: T_2	11
5.2 Beeswax: T_2	12
6 Feb 26, 2026	13
6.1 High-strength Neoprene: T_2	13
7 March 3, 2026	14
7.1 High-strength Neoprene: T_1	15
7.2 High-strength Neoprene: T'_2	16
8 March 5, 2026	17
8.1 X-ray diffraction	17
8.2 Fluorescence	18
9 March 10, 2026	18
10 March 12, 2026	20
11 March 24, 2026	20
12 March 31, 2026	21
13 April 9th, 2026	22
14 April 14, 2026	22
15 April 16, 2026	22
16 April 21, 2026	22

1 Feb 5, 2026

Link to the lab website: <https://you.stonybrook.edu/phy445/experiment-overview/>

Note: We are allowed to have 1–2 extra days for this experiment as there was no lab manual given.

1.1 Experiment: Nuclear Magnetic Resonance

Experiments that we will run:

- We will test T_1 , T_2 , T_2^* , and T_2' for three samples: distilled water, tap water, and water with metal in it (or honey if we do not have that);

Sub-experiment one: distilled water

We got the distilled water from the HEP-ex water room; poured it into the little pipe without a pipette because the only one we found was super dirty and probably full of AIDS.

Distilled water: T_1 Experimental set-up and tuning for measuring T_1 for distilled water:

- RF frequency matching: Tuned it until there was a minimal sine curve. Measured frequency is 21.01291 MHz.
- $\pi/2$ pulse: $1.12 \pm 0.005 \mu s$
- π pulse: $7.32 \pm 0.005 \mu s$

Now we tried to measure T_1 for distilled water, but this was very hard:

- We didn't know how to export the osc. data to fit it.
- We didn't know if the pulse was just one peak of the sine curve or if it was the entire curve.
- Corliss came over and helped us kind of. We didn't know how to do anything. He pointed out that we were syncing on A when we probably should have been syncing on B .
- TA came over and helped us more; got back to my original question which I'm just going to email Prof. Liu.
- After class, emailed Prof. Liu the following:

Hi Professor Liu,

My partner and I are trying to measure the T_1 timescale for the NMR experiment in PHY 445. We know that we have to fit a damped sine curve to the response of the nucleus to the A pulse to determine T_1 . After asking both the TA and Professor Corliss for advice on how to do this, we are still not sure how to take the oscilloscope data and translate it into a form that we can fit a curve to. Do you have any advice on how to do this?

Also, you previously said that some sort of iron solution is commonly used in this experiment to test different responsivities to the external magnetic field and RF pulses. Will this be made available to us next Tuesday?

Thanks,
Will Lancer

2 Feb 10, 2026

- Prof. Liu responded:

Hi Will,

Prof. Xu Du is taking care of the NMR measurement this semester. You probably confused me as Prof. Du?

Regarding the data, usually you can take the oscilloscope data from a USB drive and import it into the computer. Then you can fit it from there. (Depending on which oscilloscope you are using this semester.)

Everything should be in the test kit – usually I will suggest you use the mineral oil first, and then the rubber.

- Whoops! Literally everyone in the lab thought that Prof. Du was Prof. Liu. #rip. #duNotLiu.
- Anyway, time to continue with this experiment.

Ok, we're going to switch our experiment to mineral oil on Prof. Liu's recommendation.

2.1 Sub-experiment: mineral oil

We're using the light mineral oil. Time to re-measure all of our stuff!

- RF frequency matching: 21.11319 ± 0.00001 MHz.
- $\pi/2$ pulse: (**didn't take for some reason?**) $\pm 0.005 \mu\text{s}$.
- π pulse: $6.76 \pm 0.005 \mu\text{s}$.

We think that we want an A pulse and then a B pulse. The former being a π pulse and the latter being a $\pi/2$ pulse. But we don't really see that on the oscilloscope. We varied tons of numbers, and then figured out that we can put just one B pulse. We did that, but when we vary τ , we don't see the amplitude of B change; this isn't what we expected, as we'd expect the curve to follow the equation $\sim 1 - e^{-\tau/T_1}$.

2.2 Enter Professor Du, the King of NMR™

We called Prof. Du over and he determined after some debugging that we were way off our resonance frequency. We were using the previous value, but he determined that it should be 21.00854 ± 0.00001 MHz. The A pulse and B pulses were wrong because we were off resonance: the

new ones are $A = 6.72 \pm 0.005 \mu\text{s}$ and $B = 2.60 \pm 0.005 \mu\text{s}$.

NOTE: ALL OF THESE VALUES ARE FOR BEESWAX, NOT MINERAL OIL.
Xu “The Goat of NMR” Du switched us to that material in his infinite wisdom.

The lesson from this is that the vertical scale of the blue line was *super* zoomed out; because of this, I tuned it to what seemed tiniest but it was wayyyy off. So, make your vertical and horizontal scales fit the scales at hand or else you will get cooked.

Great. I then was curious about what the period was and why moving it around shifted the amplitude of the B pulse, so we called over Prof. Du again and he explained it to us. We then increased the period from 100ms to 600ms because the relaxation time of oil is kinda long.

After switching the period, I noticed the frequency was looking sus and a little bit beatyastic, so we changed it again to make it not beating:

- RF frequency matching: $21.00740 \pm 0.00001 \text{ MHz}$.
- π pulse: $6.76 \pm 0.005 \mu\text{s}$.
- $\pi/2$ pulse: $3.36 \pm 0.005 \mu\text{s}$.

3 Feb 12, 2026

Got into lab. Eric and friends did some bullshit to the scope and the whatever frequency setter thing. We fixed it, and switched the RF frequency to: $21.01509 \pm 0.00001 \text{ MHz}$.

We are now going to determine T_1 for the light mineral oil. See Alo’s lab notebook for the relevant table.

3.1 Beeswax T_1 measurement

Remeasured the relevant values:

- RF frequency matching: $21.01183 \pm 0.00001 \text{ MHz}$.
- π pulse:
- $\pi/2$ pulse:

3.2 Nevermind, fml

We realized that our previous A and B numbers were for BEESWAX, NOT MINERAL OIL. So we need to retake our data!!! So fun! Experiments are just so wonderful and not absolutely shit, I just love them! This definitely isn’t wasting my life!

Relevant values:

- RF: $21.01155 \pm 0.00001 \text{ MHz}$.
- π : $7.36 \pm 0.005 \mu\text{s}$.
- $\pi/2$: $3.68 \pm 0.005 \mu\text{s}$.

Asked my best buddy ChatGPT to plot the table for me. He gave us these tables:

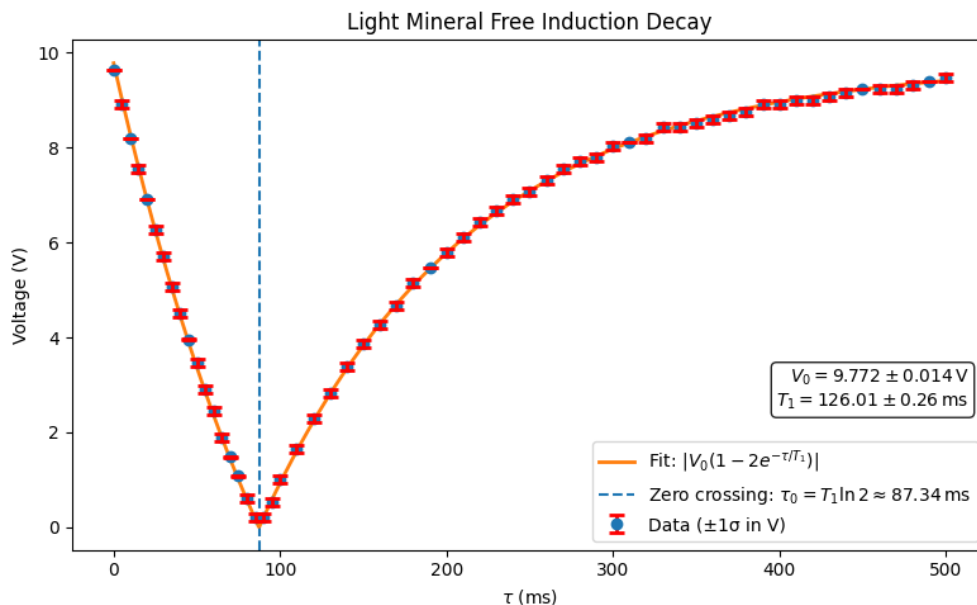


Figure 1: chungus 1

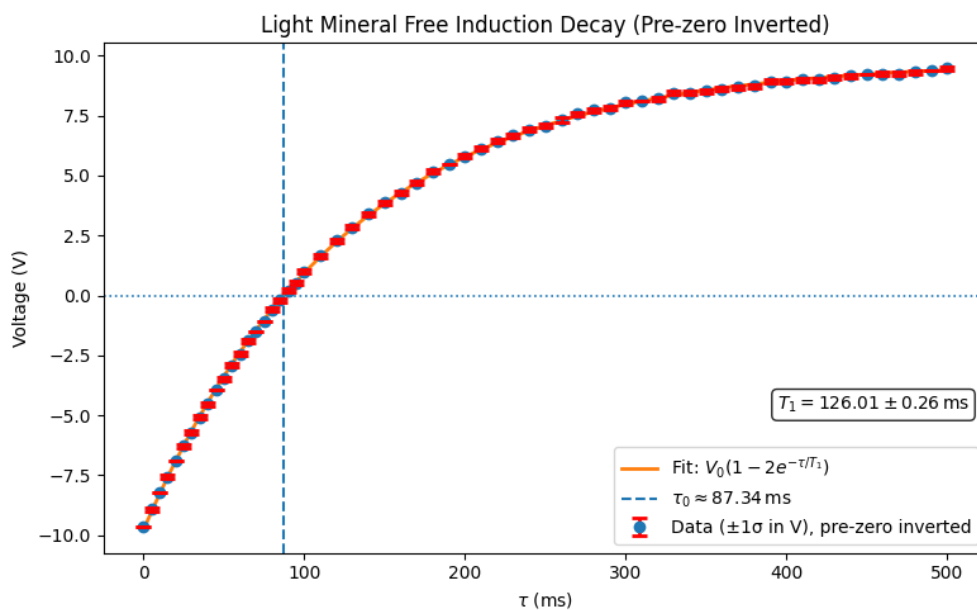


Figure 2: chungus 2

3.3 Beeswax, pt. 2

RV:

- RF: 21.00852 ± 0.00001 MHz.

- π : $6.80 \pm 0.005 \mu\text{s}$.
- $\pi/2$: $3.46 \pm 0.005 \mu\text{s}$.
- Period: $700 \pm 0.5 \text{ ms}$.

See Alo's table for the raw data.

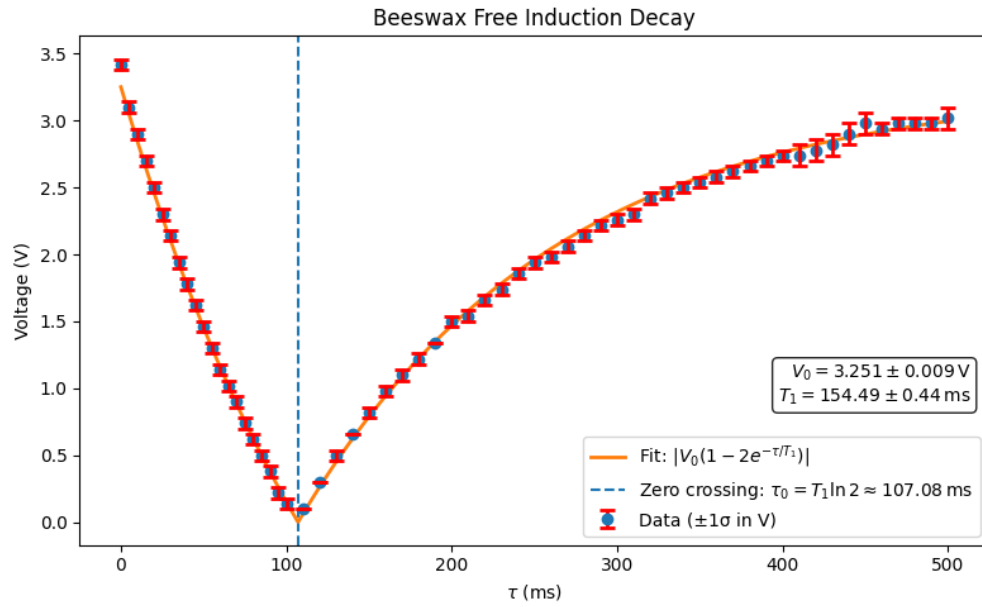


Figure 3: chungus 3

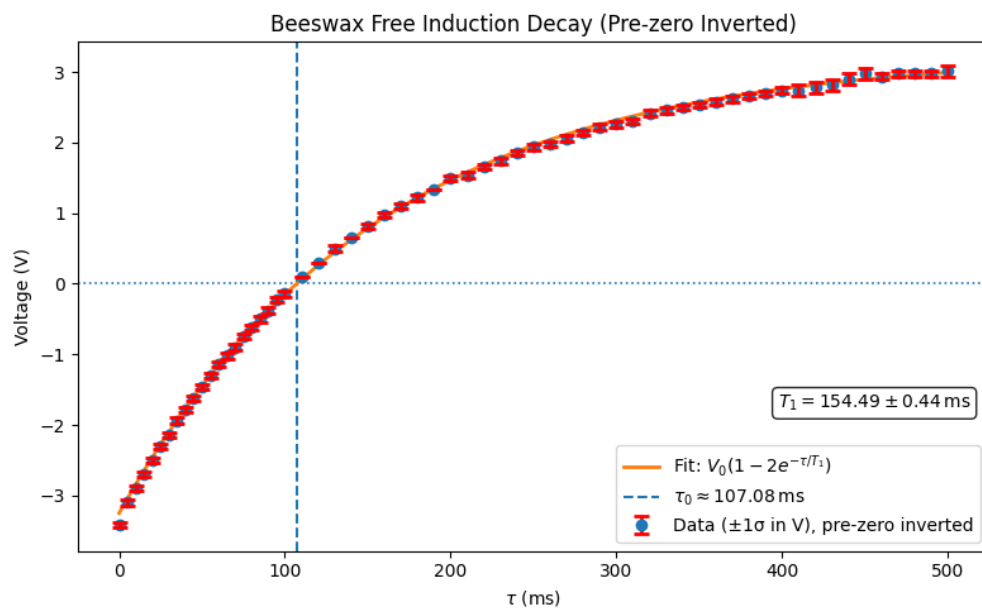


Figure 4: chungus 4

4 Feb 17, 2026

Let's chung our way through this day... this class genuinely sucks. Whatever. Let's do the spin-echo experiment.

4.1 Spin-echo: light mineral oil

- RF: 21.015224 ± 0.00001 MHz.
- π : 7.62 ± 0.005 μ s.
- $\pi/2$: 3.98 ± 0.005 μ s.
- Period: 600 ± 0.5 ms.

Señor Du helped us with these values.

Now we are going to do the spin-echo experiment. This procedure goes as follows:

1. Apply a $\pi/2$ pulse: this is the A pulse.
2. Apply a million π pulses: these are the B pulses.
3. The spin-echo curve envelope will go like e^{-t/T'_2} ; you just extract this data and fit it to a curve to get T'_2 .

We will use 30 π -pulses and fit the curve to get T'_2 . I think we will have to use thumbdrive to export the data on the scope for this one—just taking the data points may not be enough, or at least that's not what the other groups did.

We got this epic data from the spin-echo experiment:

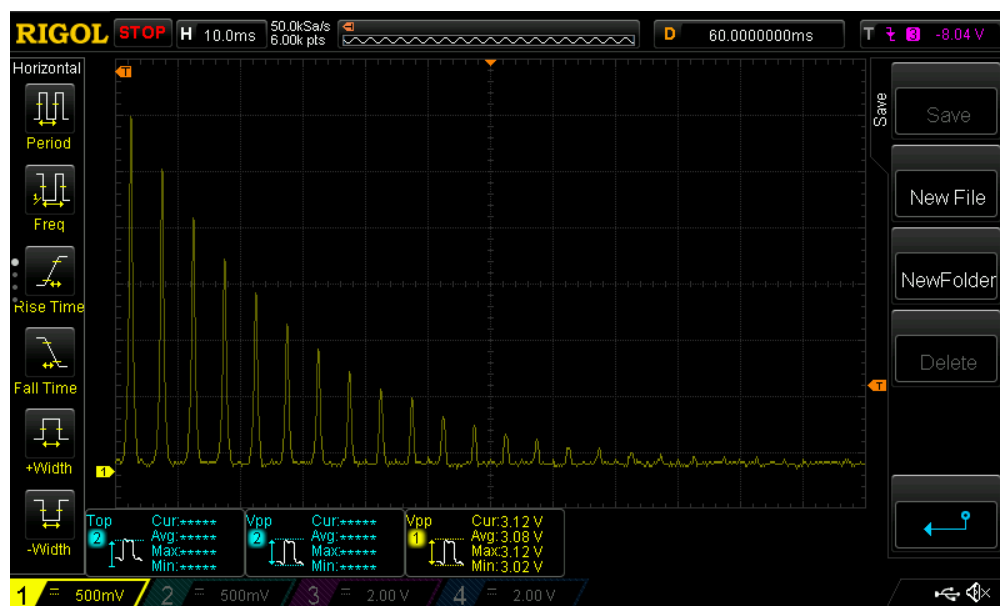


Figure 5: Spin-echo raw data for light mineral oil

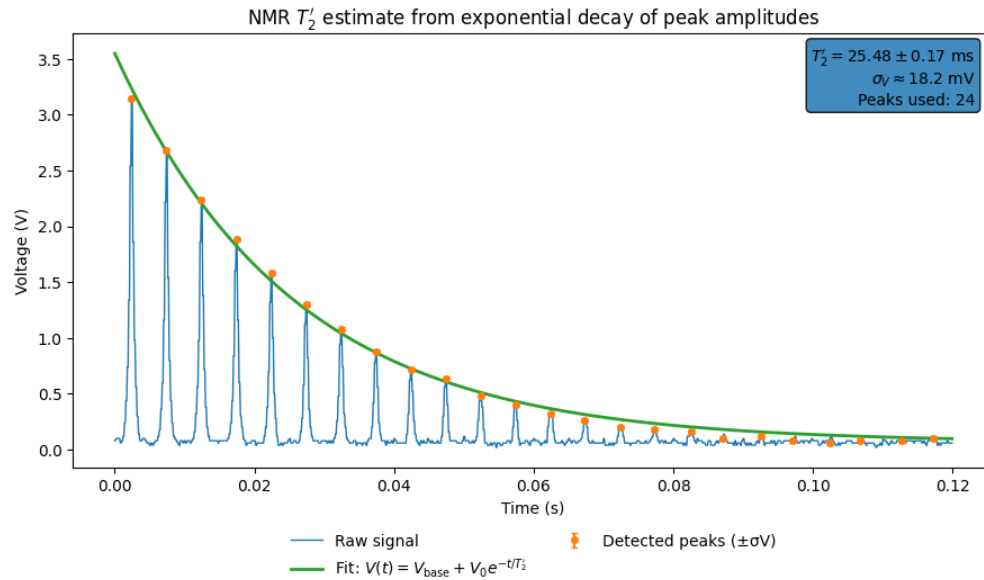


Figure 6: Spin-echo data for light mineral oil

We are the goats of NMR.

4.2 Beeswax-maxxing

Relevant values:

- RF: $21.01583 \pm 0.00001 \text{ MHz}$.
- $\pi/2$: $7.62 \pm 0.005 \mu\text{s}$.
- π : $3.98 \pm 0.005 \mu\text{s}$.
- Period: $700 \pm 0.5 \text{ ms}$.

Tables type shit!

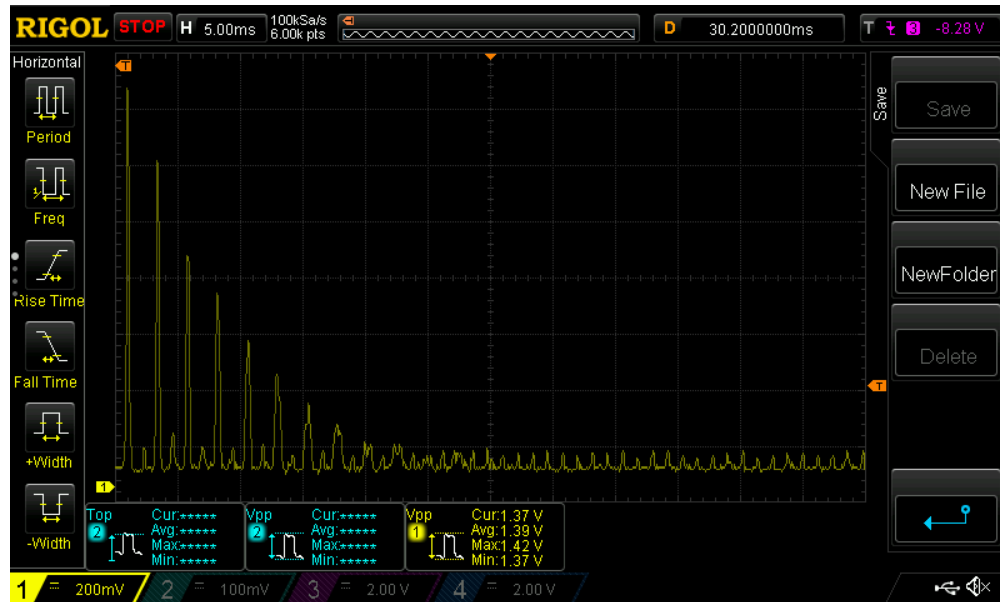


Figure 7: Spin-echo raw data for beeswax

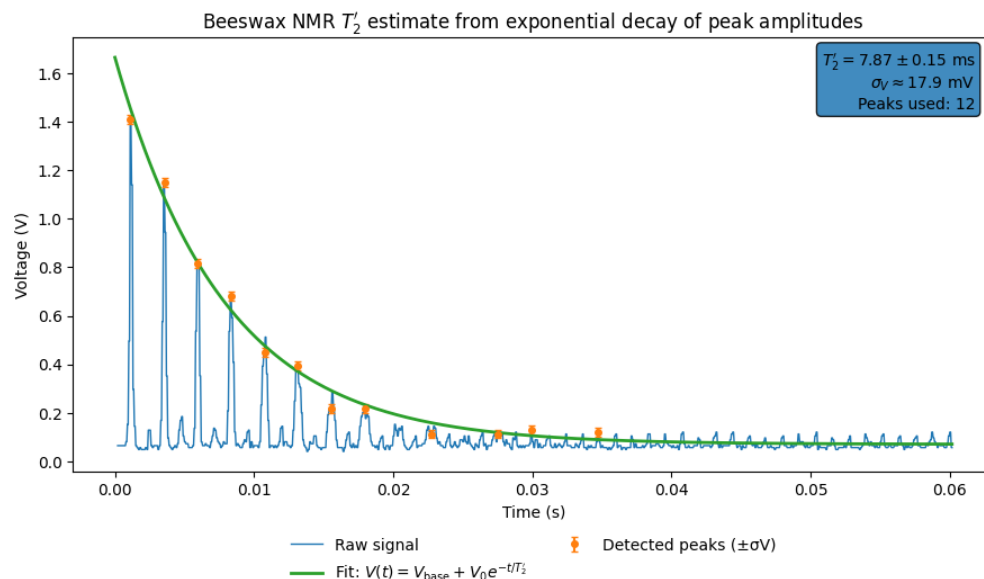


Figure 8: Spin-echo data for beeswax

We are the goats of NMR, part 2.

5 Feb 19, 2026

Back at it again. We just need to measure T_2 today, and then we are free from the shackles of NMR. The midlab report is due today; Alo submitted it for us. She now knows how to use Git, which is great.

I believe the way you measure T_2 is by pulsing a $\pi/2$ pulse and just waiting for the voltage to decay back to zero. Then you fit that to e^{-t/T_2} and you're done. Let's check with Du to see if that's Tru.

We checked with Du: this is indeed Tru. Time to Mu-ve on!

5.1 Light mineral oil: T_2

The ppl we share this experiment with messed up the stuff on the scope and with the light mineral oil sample. We now have to clean up their mess.

Relevant values:

- RF: 21.01596 ± 0.00001 MHz.
- $\pi/2$: 4.08 ± 0.005 μ s.
- Period: 600

We had to put the ring on the sample quite low to avoid a big shoulder. This wasn't the case before, but the other people on this experiment completely messed up the same and didn't store it correctly, so the light mineral oil wasn't concentrated entirely at the bottom of the tube. This lowered the measured voltage, but shouldn't be a big deal.

We got this picture from the scope, and this curve was fit to it by my bestie ChatGPT 5.2 Extended Thinking. Samriddhi was goats and let us use her personal flashdrive to get this data, and the previous data.

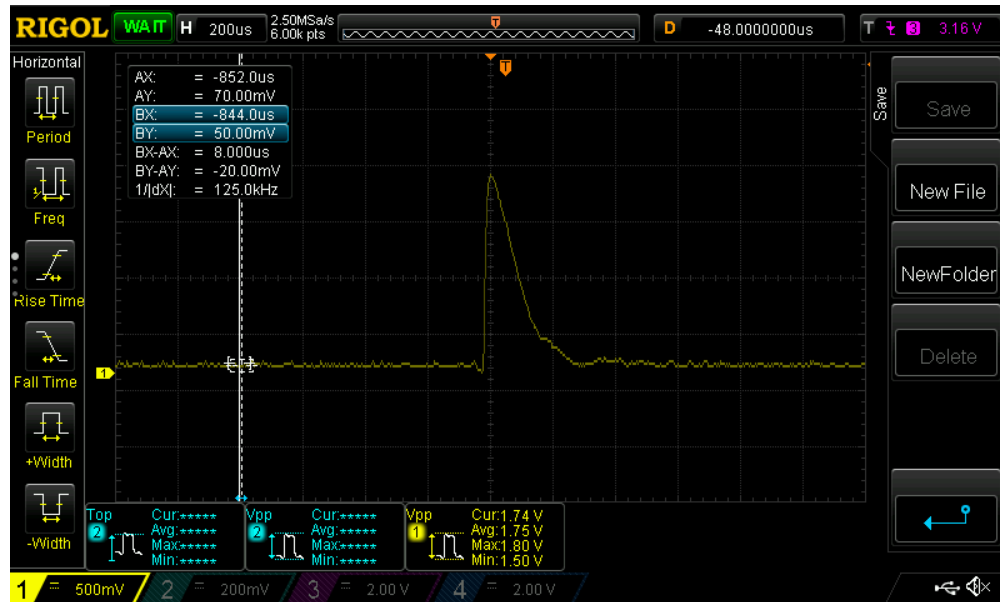


Figure 9: Raw T_2 data for light mineral oil

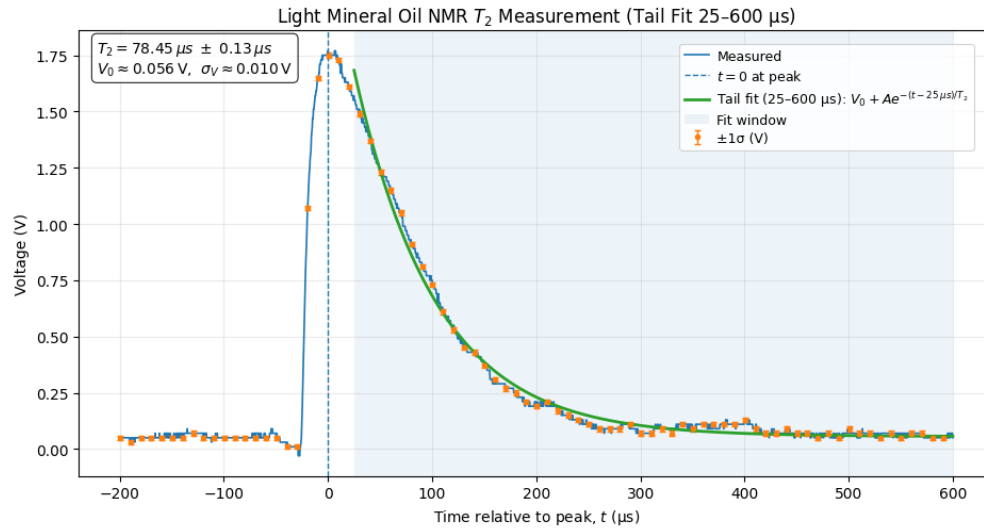


Figure 10: T_2 data for light mineral oil, fit to e^{-t/T_2} curve. $T_2 = 2.15 \pm 0.05$ ms.

5.2 Beeswax: T_2

Relevant values:

- RF: 21.01245 ± 0.00001 MHz.
- $\pi/2$: 4.10 ± 0.005 μ s.
- Period: 700 ± 0.5 ms.

Note that each horizontal block on the raw data below is 50 microseconds long.

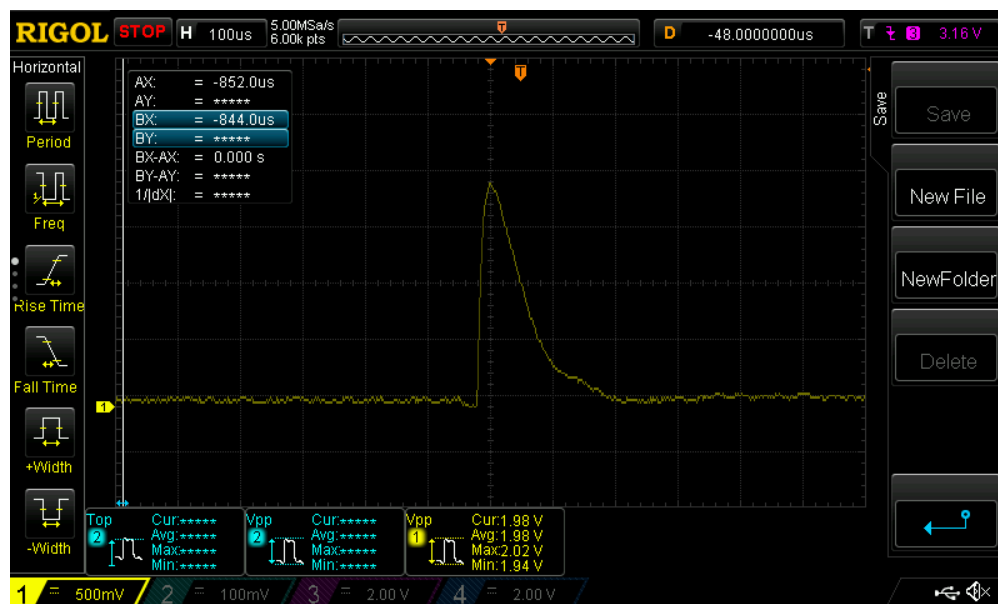


Figure 11: Raw T_2 data for beeswax

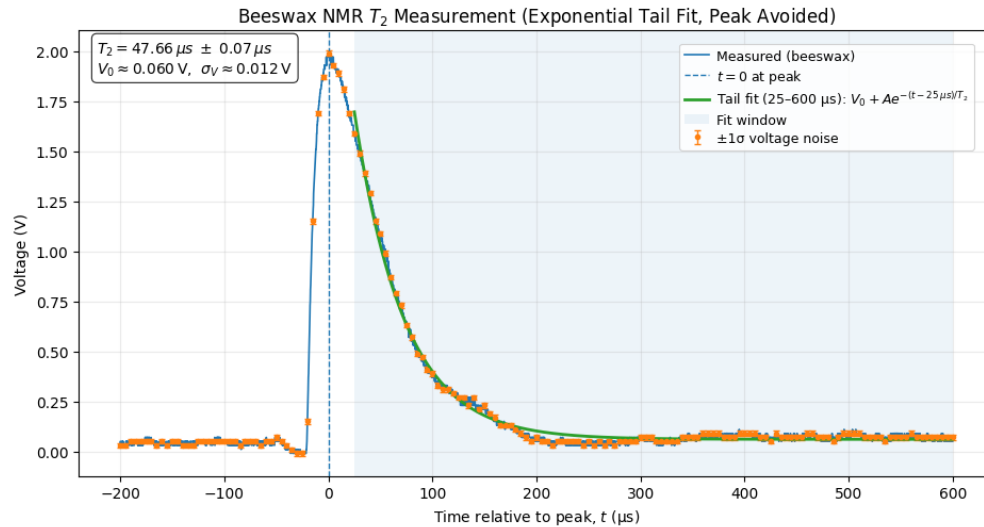


Figure 12: T_2 data for beeswax, fit to e^{-t/T_2} curve

6 Feb 26, 2026

We are back at it again. We are done with the experiment. Du encouraged us to do more samples, so we are going to do that now.

We are going to do high-strength Neoprene (129N13).



6.1 High-strength Neoprene: T_2

We are going to do the T_2 measurement for high-strength Neoprene now. Recall that the way to do this it to apply a $\pi/2$ pulse and export the graph and fit it to e^{-t/T_2} .

We took the picture but then I pressed a button that messed up literally everything on the scope. So had to call Du over and retake things.

Relevant values:

- RF frequency matching: 21.01270 ± 0.00001 MHz.
- $\pi/2$ pulse: 3.40 ± 0.005 μ s.
- Period: 650 ± 0.5 ms.

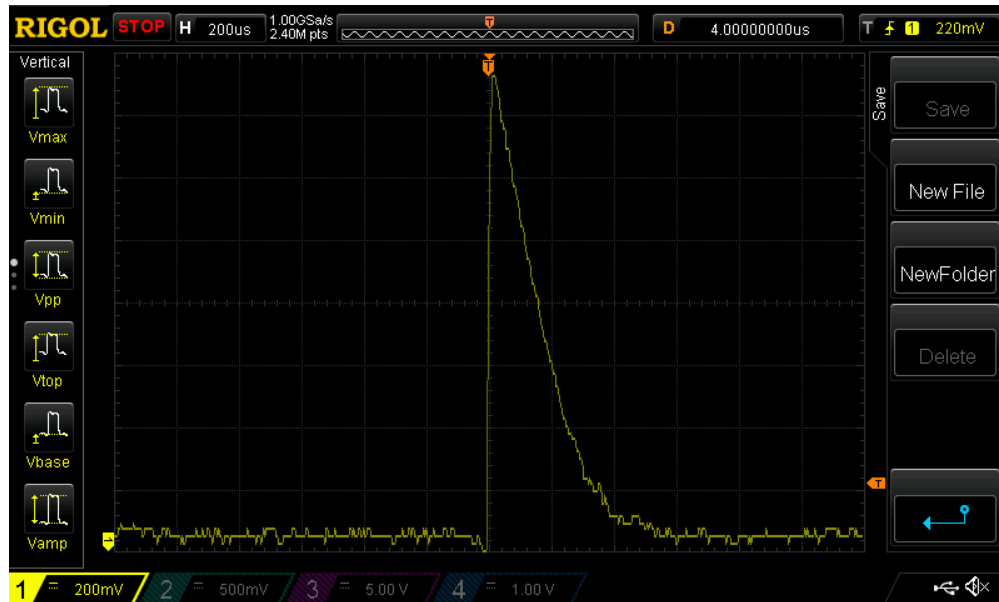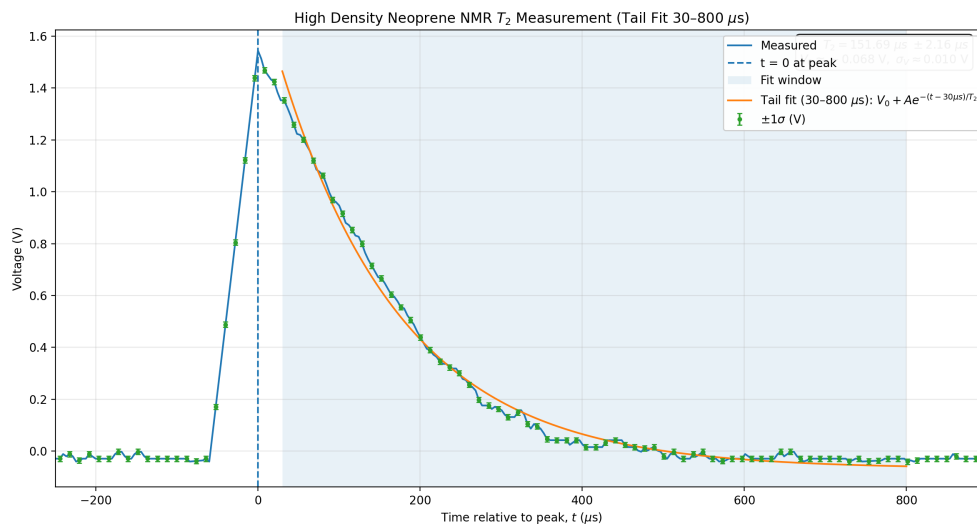Figure 13: Raw T_2 data for high-strength NeopreneFigure 14: T_2 data for high-strength Neoprene, fit to e^{-t/T_2} curve.

We tried to measure T_1 for the high-strength Neoprene, but the scope's trigger was super weird. We tried to fix it for a while, and both of the TAs did as well. It didn't work. We tried to get Du, but he was gone. So we just left.

7 March 3, 2026

WE BACK AT IT AGAIN!!! Du fixed the scope for us in about two minutes, as always (he's the goat), and so now we're going to measure T_1 and T_2' for the high-strength Neoprene.

7.1 High-strength Neoprene: T_1

Relevant values:

- RF frequency matching: 21.00766 ± 0.00001 MHz.
- $\pi/2 = B$ pulse: 3.38 ± 0.005 μ s.
- $\pi = A$ pulse: 6.74 ± 0.005 μ s.
- Period: 650 ± 0.5 ms.

Recall the procedure for T_1 measurement:

- Find A and B pulses etc.
- Find total τ for it to return to its max value.
- Take ~ 50 data points, measuring the voltage at each time
- Fit the data to the $\sim (1 - e^{-\tau/T_1})$ curve; invert the data pre-zero to get a nice damped exponential fit.

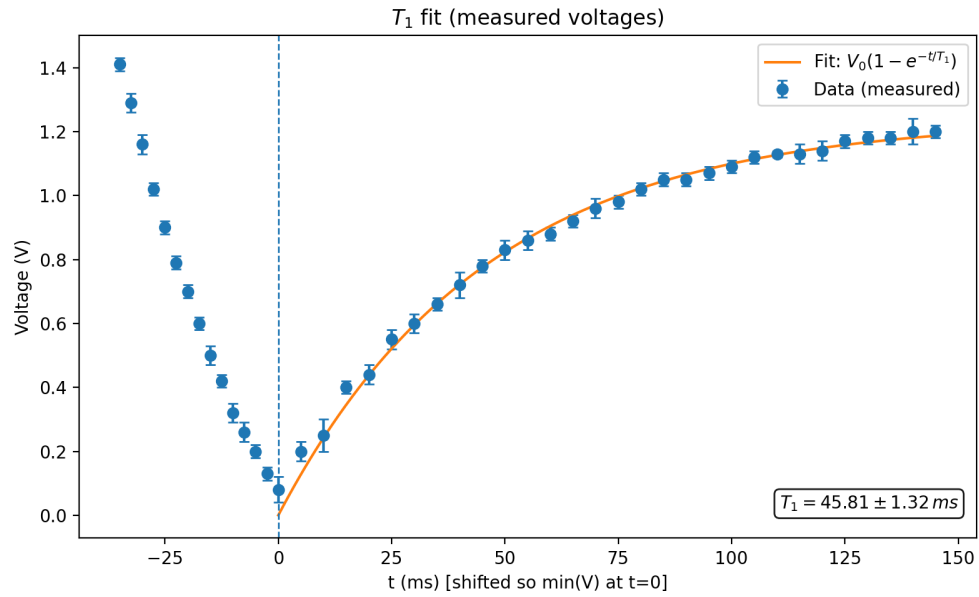
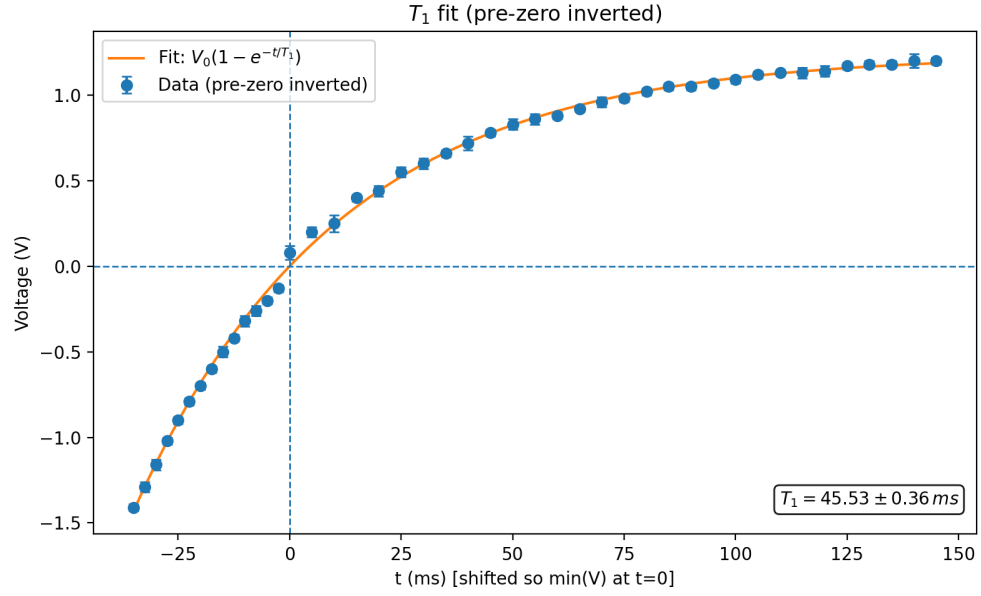


Figure 15: T_1 curve with no inversion



Figure 16: T_1 curve with inversion

7.2 High-strength Neoprene: T'_2

Relevant values:

- RF frequency matching: $21.00766 \pm 0.00001 \text{ MHz}$.
- $\pi/2 = A$ pulse: $6.74 \pm 0.005 \mu\text{s}$.
- $\pi = B$ pulse: $3.38 \pm 0.005 \mu\text{s}$.
- Period: $650 \pm 0.5 \text{ ms}$.

Did the normal fits.

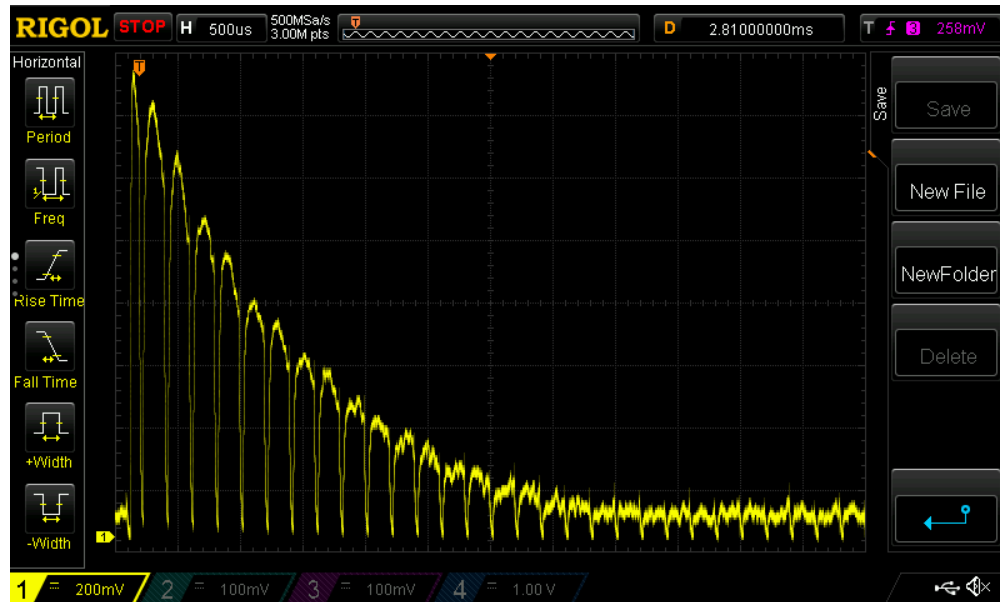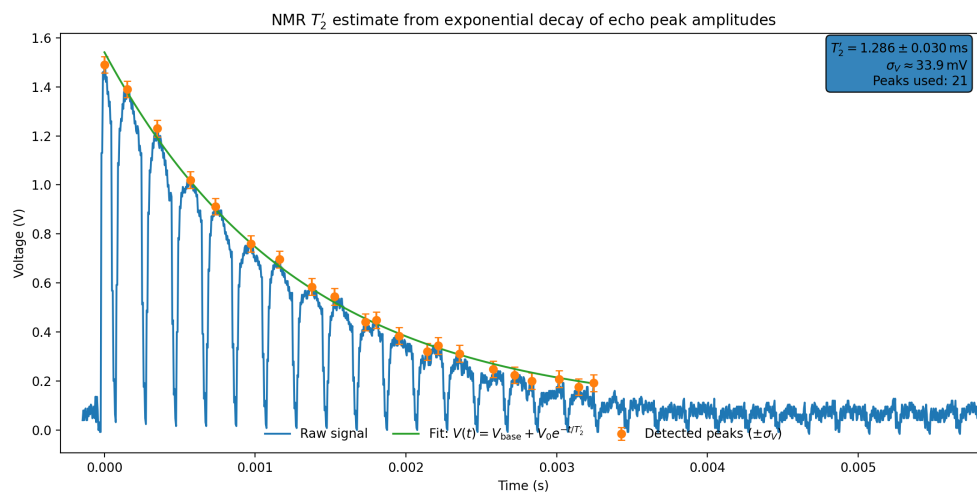Figure 17: T_2' neoprene raw dataFigure 18: T_2' neoprene fitted data

8 March 5, 2026

We no longer have the Du by our side... please take a moment of silence. We have Corliss, though, and he's great, so we're chilling.

We're starting on X-ray diffraction and fluorescence today. We had a discussion with Corliss to start the lab, and then Kevin set us up with the equipment.

8.1 X-ray diffraction

Dis lab easy af!!!

We swept the following voltage values:

- 35 kV
- 34 kV
- ...
- 24 kV

in one kV steps.

8.2 Fluorescence

We replaced the detector with the other detector. Then adjusted the detector to be at 90° to the salt crystal, and measured the following distances between the emitter-sample and sample-detector:

- Emitter-sample: 5.5 cm
- Sample-detector: 5.5 cm

We then adjusted the salt crystal from before to be at 45° to the detector.

We replaced the salt crystal with the FeZn sample. We did the set-up detailed in the lab manual and collected data for the FeZn sample. We will do the rest of fluorescence next time.

9 March 10, 2026

I'm sick today, so Alo is doing all of the equipment-touching. I'm just coaching her through it and doing the data analysis.

Alo set up the X-ray fluorescence set-up, adjusted things to be 5.5cm, and then started things up.

I had Chat 5.4 Extended Thinking fit the data to the set of linear curves we want and make a table from the data as well. We got

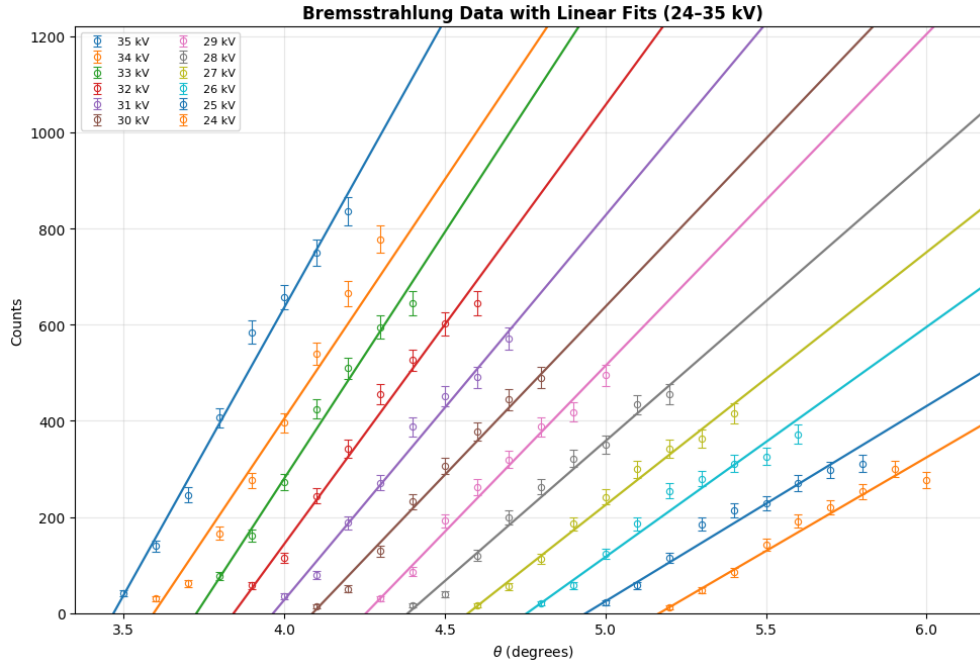


Figure 19: Linear fit lines for 35 – 24kV sweeps.

Voltage (kV)	Slope (Counts/deg)	σ_{Slope}	y-intercept (counts)	$\sigma_{\text{intercept}}$
35	1196.95	25.66	-4152.67	93.10
34	995.83	23.03	-3579.30	85.65
33	1025.46	31.43	-3822.32	124.35
32	913.26	24.57	-3509.83	99.90
31	799.99	21.64	-3172.94	89.61
30	699.66	18.42	-2860.99	77.26
29	690.37	20.15	-2937.61	89.51
28	582.29	15.26	-2553.74	69.07
27	526.24	14.66	-2406.54	69.39
26	478.72	14.58	-2276.89	72.22
25	406.30	13.93	-2007.01	72.06
24	389.65	12.63	-2013.80	67.37

Table 1: Fitted Parameters and Uncertainties for Each Voltage

We took fluorescence data, but then realized that it was around half the size that it should be (i.e. half the number of counts), but its shape is qualitative what we'd expect (four peaks; two sets of two peaks, one big, one small).

We asked Kevin O'Shea why this was, and he said it may be because we are using the collimated emitter. We switched out our collimated emitter for the circular one, and re-ran our data on the FeZn sample.

It is much better. Its counts are very high, around 3x of what they were before.

10 March 12, 2026

We retook FeZn. Fit its values to determine its peaks. Used Chat to do so—I ran out of file uploads so Alo did it for us.

Then we took the following samples:

- Ti — titanium
- Ag — silver
- Rh — rhodium
- Au — gold

11 March 24, 2026

Long time no see—it was nice while it lasted.

Chat analyzed our data from last time and produced the following plots:

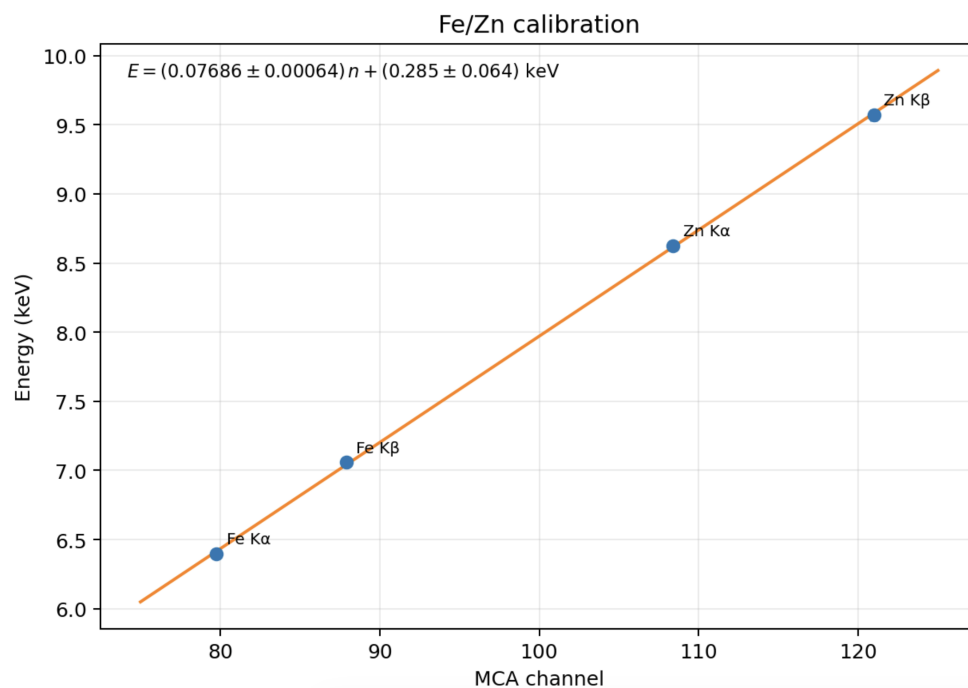


Figure 20: Fluorescence calibration data for FeZn

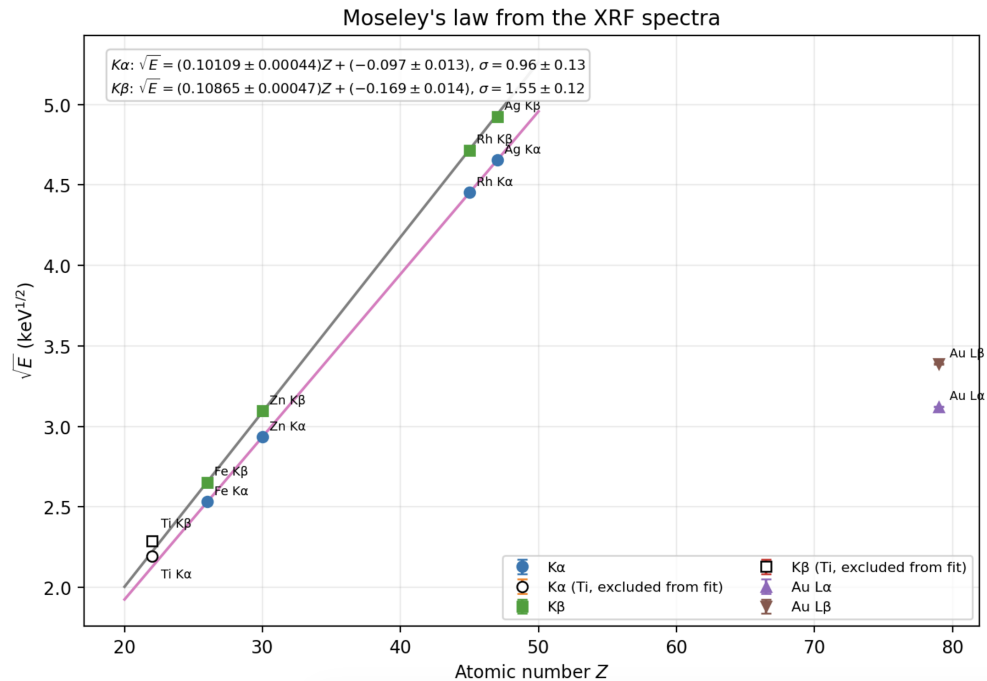


Figure 21: Fluorescence data for other samples

12 March 31, 2026

Corliss told us to retake our fluor data for gold and to take more samples to see if the trend of the L_α and L_β lines only appearing for higher atomic numbers holds. So we measured the following samples, measuring the reference FeZn again.

- Gold
- Lead
- Uranium
- Tungsten
- FeZn

These are ordered by atomic number as

- FeZn
- Tungsten
- Gold
- Lead
- Uranium

See the analyzed data in the assets folder: it's the only PDF in there.

13 April 9th, 2026

Starting Mossbauer today. Waited around an hour for Prof. Drees to come help us get set up. Once he came over, we learned a bit about the experiment and got the basic things set up. See the pictures on your phone for the exact set-up.

We then proceeded to find PPE in the form of gloves that we got from the laser labs downstairs. We then did the first part of the experiment, i.e. the stuff on page 2 of the detailed instructions (as in, after the page count resets at page 11—page 2 of that part).

Specific comment:

- The left triangle should be around 227 and the right triangle around 350. You can see the csv files to recalibrate this on tuesday. We do not have enough time to use these today.

14 April 14, 2026

We did the calibration for the iron sample yesterday, and today we are going to do some of the other experiments. We are waiting for Drees to bring the experiments, as we don't know if we have to run them overnight or something.

We set the first experiment running and waited around thirty minutes to verify the correct behavior. It looked good so we left.

15 April 16, 2026

Just started the data and left.

16 April 21, 2026

Did all of the data analysis today for our previous experiments. Didn't start any new experiments because we didn't know how + Drees wasn't here.

We did the Fe57 analysis, the enriched iron analysis, and the isomer shift analysis today.